




# 15 & 50mL Polypropylene Tubes

Disposable polypropylene tubes and caps with recessed lid and sealing ring to prevent leaking. Frosted writing surface allows for permanent identification and labeling. Autoclavable to 121°C and can be stored in freezers to -80°C. Manufactured in ISO 9001:2000 and 100,000 grade cleanrooms. Sterile.

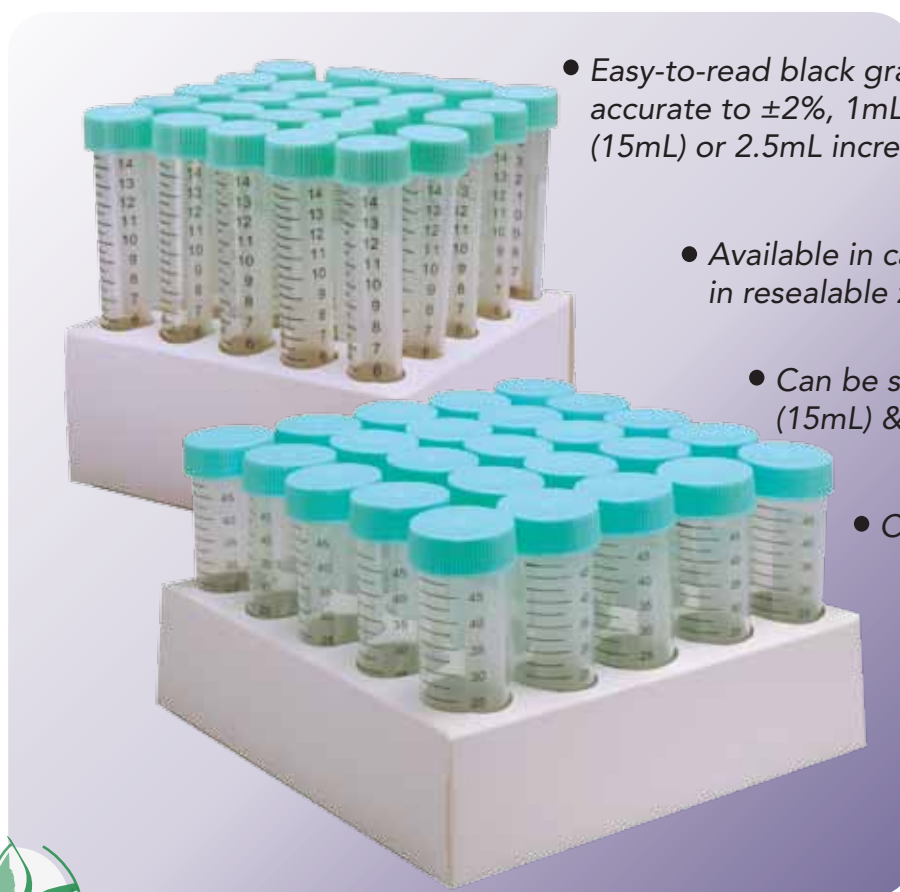

- Easy-to-read black graduations are accurate to  $\pm 2\%$ , 1mL increments (15mL) or 2.5mL increments (50mL)


- Available in cardboard racks or in resealable zip closure bags

- Can be spun up to 12,000g (15mL) & 15,000g (50mL)

- Conical Bottom

- Sterilized with gamma irradiation, RNase & DNase Free, DNA free, Non-pyrogenic

- Chloroform resistant graduations and writing area

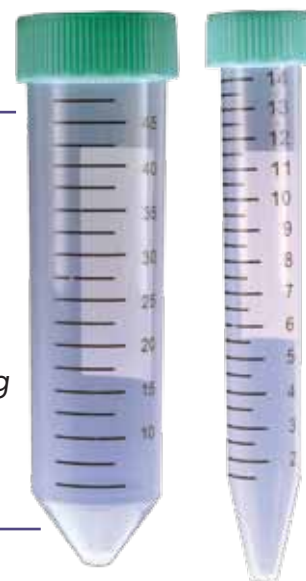

Environmentally friendly, fully recyclable cardboard racks

Catalog No.	Capacity	Packaging	Qty.
TC1500	15mL	Resealable Bag	25/Bag, 500/Case
TC1500R	15mL	Cardboard Rack	25/Rack, 500/Case
TC5000	50mL	Resealable Bag	25/Bag, 500/Case
TC5000R	50mL	Cardboard Rack	25/Rack, 500/Case



Ph: 800-886-8675 • Fx: 847-622-0454

[www.argos-tech.com](http://www.argos-tech.com)

Distributed by: