



## Disposable Plastic Serological Pipettes

Designed for accurate and rapid fluid dispensing; these polystyrene pipettes are sterile, pyrogen free and utilize a unique filter plug to prevent overfilling. The filter plug eliminates lint fibers and ensures plug stays in place. The sharp, black, highly legible permanent graduations prevent ambiguity and guarantees precise reading and accurate dispensing to +1 - 2%. Reverse graduations show volume removed or remaining on all pipettes larger than 1mL. Generous negative graduations provide additional pipetting volume.

Argos pipettes are made from transparent medical grade polystyrene sterilized by gamma irradiation. Each lot is tested for endotoxins (C.5 Eu/mL) (ISO T0993-11), cytotoxicity (ISO10993-5) and are guaranteed non-hemolytic (ISO10993-4). Ideal for tissue culture applications.

Pipettes are supplied individually wrapped in a fiber-free paper & plastic combination that is easy to open or available in bulk. The uniform mouth piece fits all major pipettors and have color coded TD rings for ease of identification. 1.1 & 2.2 bacteriological pipettes meet the A.P.H.A. standards for the examination of dairy products. Calibrated TD with blow out. Meets AST < E934.


Dispenser box allows user to remove pipettes without opening top. Universal stacking design requires less laboratory space and each box is color coded for easy recognition.



Color coded rings for easy identification



Fiber free paper wrap ensures quick peel apart opening

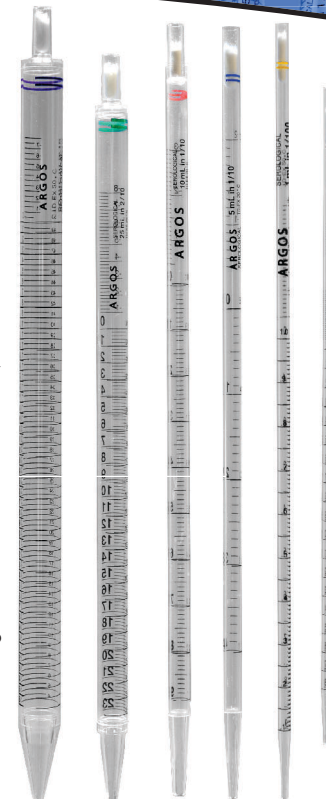

Cat. No.	Capacity mL	Graduation mL	Negative Grad. mL	Color Code	Qty
<b>Bulk Pipettes</b>					
PMB11	1.1	0.5, 1.0, 1.1	-	Purple	25/Bag, 1000/Case
PB01	1	0.01	0.01	Yellow	25/Bag, 1000/Case
PB1	1	0.1	0.1	Red	25/Bag, 1000/Case
PB2	2	0.01	0.01	White	25/Bag, 1000/Case
PB5	5	0.1	0.1	Blue	25/Bag, 500/Case
PB10	10	0.1	0.1	Orange	25/Bag, 400/Case
PB25	25	0.2	0.2	Green	10/Bag, 150 Case
PB50	50	0.5	0.5	Purple	10/Bag, 100/Case
PB10W	10 (Wide Tip)	0.1	0.1	Orange	25/Bag, 400/Case
<b>Individually Wrapped</b>					
PMS11	1.1	0.5, 1.0, 1.1	-	Purple	500/Case
PS01	1	0.01	0.01	Yellow	500/Case
PS1	1	0.1	0.1	Red	500/Case
PS2	2	0.01	0.01	White	500/Case
PS5	5	0.1	0.1	Blue	200/Case
PS10	10	0.1	0.1	Orange	200/Case
PS25	25	0.2	0.2	Green	150/Case
PS50	50	0.5	0.5	Purple	100/Case



Ph: 800-886-8675 • Fx: 847-622-0454

[www.argos-tech.com](http://www.argos-tech.com)